



FKG
swiss endo



NiTi rotary files
for mechanised
glide path preparation

ScoutRace

Sequence designed for canals with severe curvature

► The FKG Race range

An internationally-recognised reference, Race (Reamer with Alternating Cutting Edges) is the most complete selection of NiTi rotary files for root canal treatment.

ScoutRace

The range designed for canals with severe curvature

ScoutRace is a sequence of three instruments with .02 taper and with ISO sizes of 10, 15 and 20. They are designed to mechanically prepare root canals with severe curvature, or «S» type canals. The files operate at a speed of 600-800 rotations per minute, and are used after determining the working length. The ScoutRace range enables the glide path shaping faster and safer than

with manual instruments. The anatomy of the median and apical thirds are well respected, to experience no transportation and nor ledges. The files adapt perfectly to the canal's shape, thanks to their .02 taper and remarkable flexibility. The ScoutRace can be used without pressure up to working length. Root canal treatment can then be carried out with iRace or BioRace sequences.

Illustration of a case with a severe apical curvature

Courtesy of Dr. Gilberto Debelian (Norway)



Pre-operation
Diagnosis: pulpitis

Difficulty:
Double curve

Post-operation
Preparation: 35/.04


FKG Race: Specifications

Several exclusive technological advances lie behind Race range's efficacy and worldwide renown:

- ▶ Exclusive rounded safety tips for precise guiding
- ▶ Alternating cutting edges to eliminate screwing-in effect
- ▶ Sharp edges for best cutting efficiency
- ▶ Electro-chemical polishing for better resistance to torsion and fatigue
- ▶ SafetyMemoDisc (SMD) to master metal fatigue

ScoutRace availability

Individual Files	5 pieces/box
Assortment	2 × 10/.02, 2 × 15/.02, 2 × 20/.02
Lengths	21 mm, 25 mm, 31 mm

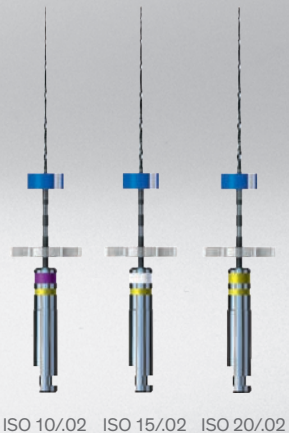



FKG Dentaire SA

Founded in 1931, FKG Dentaire SA is today at the very cutting edge of development, production and distribution of dental products for dentists, endodontists and laboratories. Precision has always been central to FKG, and the company naturally made its base in «Watch Valley», home to most Swiss watchmaking and microtechnical companies.

FKG was given a new lease of life in 1994, when Jean-Claude Rouiller took up the reins. A CEO with a vision, he set out a strategy based on innovative products and the concept of machinery specially manufactured for the dental industry. Mr Rouiller also enlarged the

distribution network to encompass more than 80 countries worldwide. FKG products are certified according to international norms and standards.

The Swiss Venture Club named FKG «Western Switzerland Company of the Year 2012», a reward for the company's dynamism, high product quality, and its continuing innovation.

**« Western
Switzerland
Company of the
Year 2012 »**



FKG Dentaire SA, Crêt-du-Loclé 4, CH-2304 La Chaux-de-Fonds, Switzerland
T +41 32 924 22 44, F +41 32 924 22 55, info@fkg.ch, www.fkg.ch